

Warning:

To reduce fire, burns, electric shock, or other injury risk, please observe the followingInstructions.

1)When nobody cares it, please pull out the plug from the socket. In addition, when it is out of using or opening it, you must pull out the plug.

2)To avoid overheat, electric shock, or injury, don't use the product in pillows or blankets under operation.

3)If the product is close to the children, physical impairment and disability when using , the need comprehensive the guide.

4)Please use it according to the manual of the product. Don't use non-manufacturing accessories recommended by other manufacturer.

Danger:

In order to reduce the risk of electric shockdanger, please follow the instructions below

1) After using or before cleaning products, please pull out the plug.

2)If the product falls into the water, don't touch the wet product. You must pull out the plug.

3)Do not use the product when taking a shower or a bath.

4)Do not place the product in water or other liquid.

Troubleshooting


If you encounter problems when using it,please refer to the following tips and solutions .If the massager is still in the fault state or not available,please turn off the power and unplug the plug.Users do not disassemble or repair massager by yourself.You should ask the seller to help you.

Condition	Possible reasons/Solutions
Unable to start the massager	Make sure the power plug is in socket. Confirm the switch is on.Confirm if the button has been pressed.
When using,it suddenly stops to work.	Please check whether the power plug is pulled out or not.Confirm if the plug has a bad connection. Maybe press the button control carelessly. Automatic timer sets for 15 minutes time has passed. If the overheating protection device starts, please turn off the main power supply,and wait for the massage cooling. If the problem persists,please do not continue to use and contact the company repair service center.

If you have any questions or problems during use, please send email to Z2685934551@hotmail.com with Amazon order number. We reply email everyday.

- ★ If a product is something wrong ,do not privately disassembly it. must pay attention to it.
- ★ The product is with heating on the surface,the people who are insensitive it will be damaged.
- ★ Please do not place the product on the solid surfaces to use ,because Do not use it in the bathroom and humid environments.
- ★ will cause the product to stop work.
- ★ Don't pull the connection by force,the bad connection of transformer
- ★ Make sure that products and transformer correctly contact.
- ★ the rated voltage and then use.
- ★ It must be provided in the transformer (DC12V 2A)or confirm

Notice

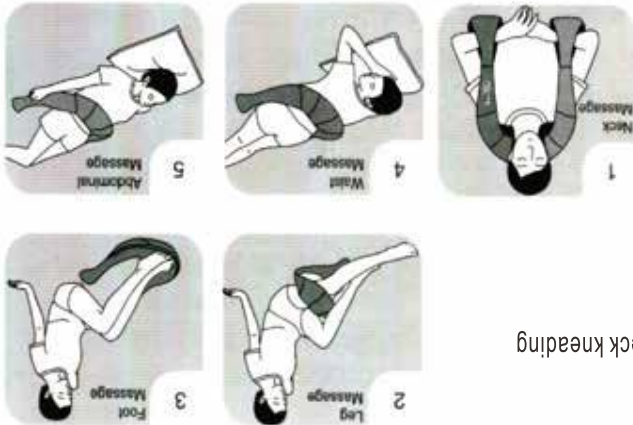

Product Name:Massager of neck kneading

Rated voltage:DC12V

Rated power:24W

Rated frequency:50HZ

Device type: ☐

Product weight: 3 Pounds

you can enjoy massage fun in any place. Equipped with external adapter,and you can use it by 12V vehicle power supply,so

04 Home used and car used

Warm function will bring you different feelings;It will bring you the warm feeling while it is having effective massage for you.

03 The unique warm function

heads knead and massage upward or downward,let you feel relaxed. Give up the traditional massage,and have a deep massage on the neck,4 massage

02 Unique way of kneading massage

The new design,exquisite curve design,tightly massage every position.

01 Massage of Neck and shoulder

Parts name and function

- HT(heater) If you want to adjust the buttons while working,press one time,it will be closed. Press again,it will be open.
- ON/OFF Working time:15mins It will be closed after using for 15mins
- DR(DIRECTION) You can adjust the direction of turning. It can be used effectively inspecific areas. You can massage specific part of your body.
- Speed Adjustment You can use it to speed up or slow down,the kneading speed:medium, fast and slow



REGULAR MASSAGE CAN IMPROVE THE WHOLE BODY FUNCTION.

- Improve the quality of sleep. Massage can make the body relax,eliminate fatigue,and relieve pressure. It also can improve sleep,and improve sleep quality.
- Improve body function Massage can promote blood circulation for the body to provide adequate nutrients,improve lymph detoxification. Relax the muscles to maintain body balance of yin and yang.
- Condition gastrointestinal function Massage can improve gastrointestinal function,promote the new supersedes the old,improve digestion. Help digestion,relieve constipation and other gastrointestinal adverse circumstances.

HT-505



We will upgrade products, if we change parameters and appearance, there will be without another notice, the product in kind prevail.

Please keep this manual,and read carefully before using.